

Instructions for Use


INTENDED USE

DNAgard[®] Saliva HT is intended for the collection and preservation of human saliva DNA. This product is For Research Use Only. Not for use in diagnostic procedures.

SUMMARY AND PRINCIPLES

DNAgard[®] Saliva HT is designed for efficient self-collection, preservation, shipping, storage, and automated processing of human saliva samples for DNA purification and analysis. Saliva DNA samples collected using DNAgard Saliva HT are stable at room temperature (15°C to 30°C) for up to 12 months, and yield high quality DNA using a wide selection of commercial DNA purification kits.

KIT CONTENTS

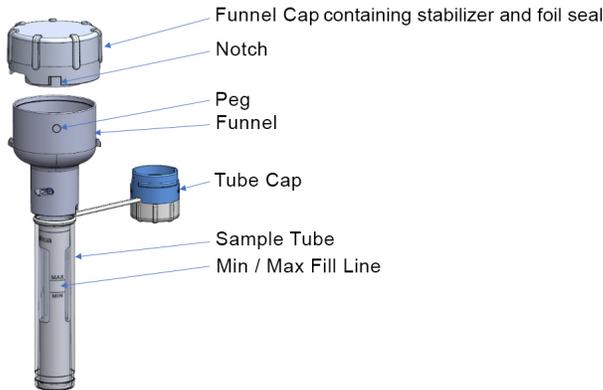

WARNINGS AND PRECAUTIONS

1. For Research Use Only. Not for use in diagnostic procedures.
2. Do not use if packaging is damaged or foil seal in funnel cap is broken or leaking.
3. Do not use beyond expiration date.
4. Avoid skin contact with DNAgard Saliva HT reagent.
5. Do not ingest DNAgard Saliva HT reagent.
6. Choking Hazard: Tube Cap is a potential choking hazard.

Contents of this reagent may cause skin and eye irritation.

1. In case of skin contact, immediately wash with soap and water and rinse thoroughly.
2. In case of eye contact, rinse immediately with plenty of water for at least 15 minutes and seek medical advice.
3. If accidental swallowing occurs, immediately seek medical advice.
4. Refer to SDS in case of accidental ingestion or skin contact. All SDS information is available at http://biomatrica.com/support_sds.php.

STORAGE AND STABILITY

Pre-Collection:

1. Store DNAgard Saliva HT kits at room temperature (15°C to 30°C) until the end of shelf life marked by the expiration date.
2. Avoid direct exposure to sunlight.


Post-Collection:

1. Store DNAgard Saliva HT sample tubes at room temperature for up to 12 months.


COLLECTION PRECAUTIONS

1. Do NOT eat, drink, smoke, or chew gum 30 minutes prior to providing saliva sample or results may be affected.
2. Do not remove the foil seal from the funnel cap.
3. If liquid does not dispense after twisting the funnel cap, press and twist further. If liquid still does not dispense, discard the intact device.
4. Do not remove the funnel cap.


DONOR COLLECTION INSTRUCTIONS

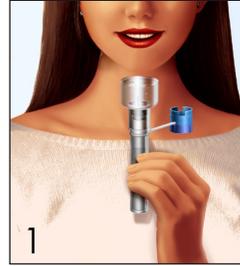


1. Hold Sample Tube upright. Spit into Funnel until saliva level (excluding bubbles/foam) is between the MIN and MAX fill lines marked on the Sample Tube. If saliva flow from Funnel to Sample Tube is slow, tap bottom of tube against a hard, level surface.

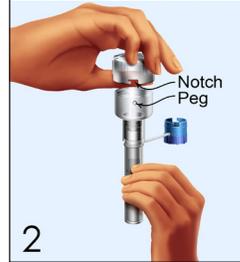


2. Place the Funnel Cap on the Funnel by aligning the Notch with the Peg. Press down firmly on the Funnel Cap to break foil seal.

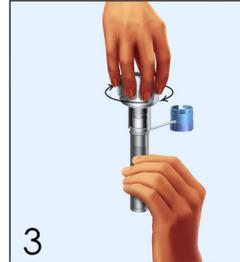


3. Twist slowly and firmly clockwise until Funnel Cap locks (clicking sound). Wait 10 seconds to ensure that all of the stabilizing solution has flowed into the Sample Tube. If flow from Funnel Cap to Sample Tube is slow, tap bottom of Sample Tube against a hard level surface.

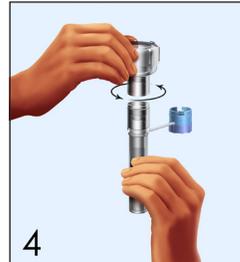

4. Unscrew and discard the Funnel and Funnel Cap together as one piece.

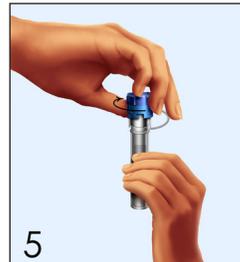

5. Screw the Tube Cap securely onto the Sample Tube.

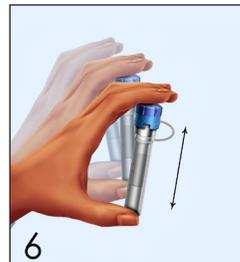

6. Shake Sample Tube for 5 seconds.

Glossary of Harmonized Symbols

REF	Catalog Number	15°C to 30°C	Temperature Limitation
LOT	Batch Code	☒	Do Not Reuse
☒	Consult Instructions	☀	Keep Out of Sunlight
☒	Manufacturer	⌚	Use By

DNAgard® Saliva *HT* Internet IFU Directions

Instructions for use can be obtained via the internet at <http://biomatrica.com/dnagard-saliva-ht.php>, by calling the following phone numbers: 1-866-379-6879 (Toll-Free in US and Canada) or 1-858-550-0308.